


Solvents Rtx[®]-1701

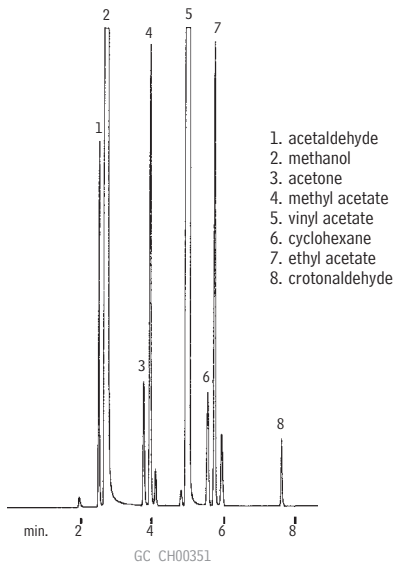

60m, 0.53mm ID, 1.5 μ m Rtx[®]-1701 (cat.# 12073)
0.4 μ L direct injection of a solvent mixture
Conc.: 2000 ppm

Oven temp.: 40°C (hold 4 min.) to 100°C @
15°C/min. (hold until last compound elutes)
Inj. & det. temp.: 250°C
Carrier gas: hydrogen
Linear velocity: 40cm/sec. (flow rate: 5cc/min.)
FID sensitivity: 16 x 10¹¹ AFS